


PASO 1: Tomar medidas

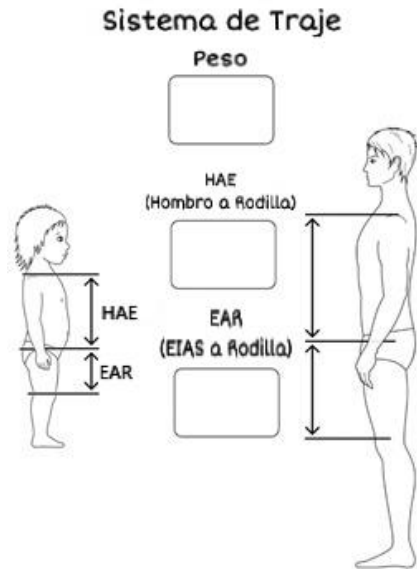


PASO 2: Talle y Selección del Sistema de Traje

Grupos de Talles	Peso (Kg)	HAE (cm)	EAR (cm)	Sistema de Traje
Infant	7 – 11	20 – 25.4	17.7 – 25.5	<input type="checkbox"/> 110
Pre-School	11 – 20,5	28 – 33	28 – 30.5	<input type="checkbox"/> 210
Pediatric	20,5 – 29,5	35.5 – 38	33 – 45.7	<input type="checkbox"/> 310
Junior	30 – 54,5	40.5 – 43	48 – 53	<input type="checkbox"/> 410
Female Med.	54,5 – 77	45.7 – 51	56 – 61	<input type="checkbox"/> 610
Female Large	77 – 91	53 – 56	63.5 – 68.5	<input type="checkbox"/> 620
Male Med.	54,5 – 77	45.7 – 48	56 – 61	<input type="checkbox"/> 710
Male Large	77 – 86	51 - 56	63.5 – 68.5	<input type="checkbox"/> 720

PASO 3: Talle y Selección de Sistema Dragonfly

Sistema Dragonfly


Largo de columna C7 a Cóccix

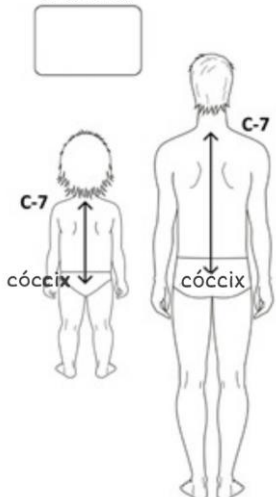

Columna (cm)	MODELO
25.5 – 30.5	<input type="checkbox"/> DFLY30
33 – 35.5	<input type="checkbox"/> DFLY33
38 – 40.5	<input type="checkbox"/> DFLY40
43 – 48	<input type="checkbox"/> DFLY48
51 - 53	<input type="checkbox"/> DFLY53

PASO 4: Medir/Seleccionar Straps de Pelvis y Tórax

Straps de Pelvis y Tórax

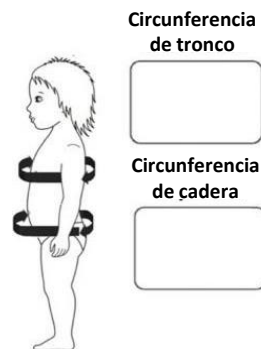

ADULTO: Strap de 5 cm

Circunferencia de Tórax (CM)	Strap de Tórax
58 - 77	<input type="checkbox"/> PS-204
78 - 102	<input type="checkbox"/> PS-206
104 - 137	<input type="checkbox"/> PS-210
139 - 157	<input type="checkbox"/> PS-214

Circunferencia de Cadera (CM)	Strap de Pelvis
43 - 59	<input type="checkbox"/> PS-204
59 - 76	<input type="checkbox"/> PS-206
78 - 104	<input type="checkbox"/> PS-210
105 - 122	<input type="checkbox"/> PS-214

PEDIATRICO: Strap de 2,5 cm

Circunferencia de Tórax (CM)	Strap de Tórax
40-61	<input type="checkbox"/> PS-103
62-71	<input type="checkbox"/> PS-105
72-84	<input type="checkbox"/> PS-107

Circunferencia de Cadera (CM)	Strap de Pelvis
38-48	<input type="checkbox"/> PS-103
50-64	<input type="checkbox"/> PS-105
66-76	<input type="checkbox"/> PS-107